

THE OXFORD® PARTIAL KNEE

The most widely used and clinically proven partial knee system in the world



References

1. Kim, K.T. *et al.* A Prospective Analysis of Oxford Phase 3 Unicompartmental Knee Arthroplasty. *Orthopedics*. 30(5 Suppl): 15-18, 2007.
2. Svard, U. and Price, A. 20-year survival & 10 year clinical results of the Oxford medial UKA. Presentation. 73rd Annual AAOS Meeting. Chicago, IL. 2006.
3. Incavo, S.J. *et al.* Tibial Plateau Coverage in Total Knee Arthroplasty. *Clinical Orthopaedics and Related Research*. 299: 81-85, 1994.
4. Goodfellow, J.W. and O'Connor, J.J. The Mechanics of the Knee and Prosthesis Design. *Journal of Bone and Joint Surgery (Br)*. 60(3): 358-69, 1978.
5. Bankston, A.B. *et al.* Comparison of Polyethylene Wear in Machined Versus Molded Polyethylene. *Clinical Orthopaedics and Related Research*. 317: 37-43, 1995.
6. Clarke, I.C. *et al.* Hip Simulator Wear Testing ArCom® vs. Extruded Bar Polyethylene. Presentation. 7th Annual Conference on Techniques & Science for Successful Joint Arthroplasty. 1995.
7. Keys, G.W. Reduced Invasive Approach for Oxford II Medial Unicompartmental Knee Replacement - A Preliminary Study. *The Knee*. 6(3): 193-96, 1999.

Surgeons are required by FDA in the United States, and strongly recommended throughout the world, to attend an Advanced Instructional Course on the Oxford® Partial Knee before using the system.

All trademarks belong to Biomet Manufacturing Corp.

This material is intended for the sole use and benefit of the Biomet sales force and physicians. It is not to be redistributed, duplicated or disclosed without the express written consent of Biomet.

For product information, including indications, contraindications, warnings, precautions and potential adverse effects, see the package insert and Biomet's website.

BIOMET®
ORTHOPEDICS, INC.

P.O. Box 587, Warsaw, IN 46581-0587 • 800.348.9500 ext. 1501
©2007 Biomet Orthopedics, Inc. • www.biomet.com
Form No. BOI0233.1 • REV113007



THE OXFORD® PARTIAL KNEE

BIOMET®
ORTHOPEDICS

KNEES • HIPS • EXTREMITIES • CEMENT AND ACCESSORIES • PMI • TECHNOLOGY

THE OXFORD® PARTIAL KNEE

- **Less** postoperative pain^{1*}
- **Less** bone removal^{1*}
- **More** natural motion^{1*}
- **Clinically Proven**
at 10, 15 and 20 years²

Tibial Component

- Anatomical shape for excellent bone coverage³



Femoral Component

- Conforming, spherical design minimizes contact stress throughout entire range of motion
- Curved inner geometry for minimal bone removal
- Multiple sizes for various patient anatomies


Mobile Meniscal Bearing

- Free-floating bearing designed to remain fully congruent with femoral component throughout entire range of motion⁴
- Only true mobile meniscal bearing knee system approved for use in the U.S.
- Increased wear resistance with ArCom® Direct Compression Molded polyethylene^{5,6}



SPECIFICALLY ENGINEERED

- Patented instrumentation** allows for easy and accurate balancing of flexion and extension gaps
- Minimal incision to avoid quadriceps disruption and allow for more rapid recovery than TKA⁷
- Size specific instrumentation allows for bone removal in precise 1mm increments and accurate, reproducible results

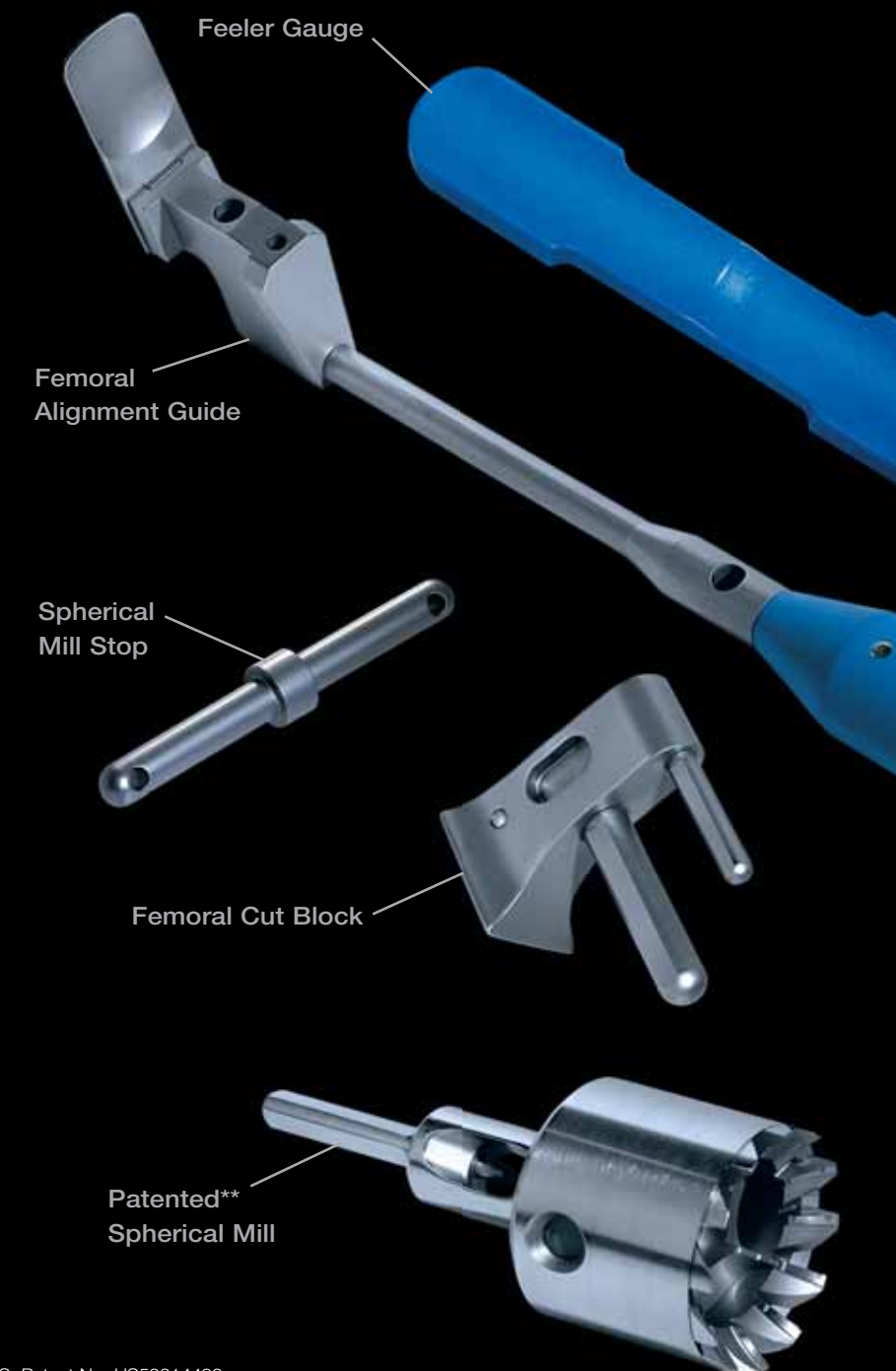

*Compared to total knee replacement

**U.S. Patent No. US53214482